



**ZEPTO™**  
PRECISION PULSE CAPSULOTOMY

---



## EXPERIENCE **ZEPTO™** **ASSISTED CATARACT** **SURGERY**

---


Precise, visually centered, and stronger  
capsulotomies in milliseconds

*The perfect capsulotomy – one less variable in refractive predictability*

**Mynosys**

# AN IDEAL CAPSULOTOMY FOR


**ZEPTO™** is a novel capsulotomy device with a disposable handpiece and tip that produces a quick, elegant, round capsulotomy in milliseconds. The tip consists of a miniature, transparent, soft silicone suction cup that houses a circular nitinol capsulotomy element, which has been refined at the micron scale to enable uniform 360-degree capsule cutting.

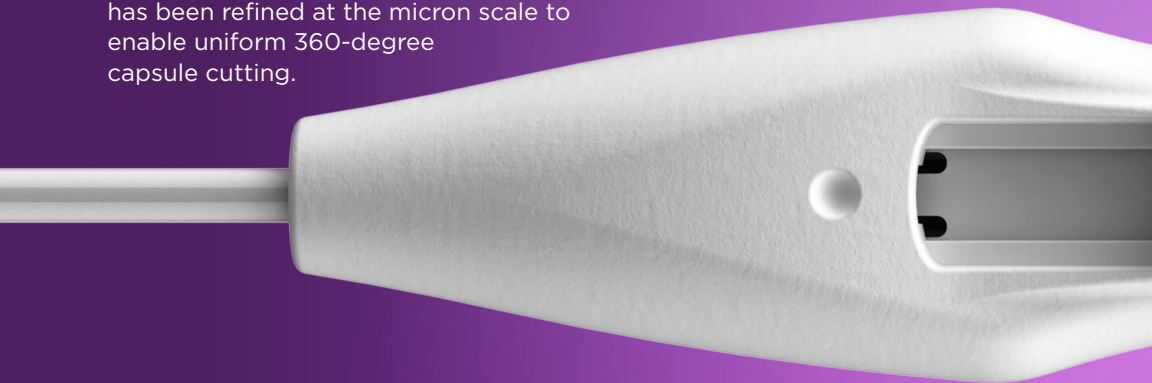


## QUICK AND EASY

Nitinol is a superelastic alloy, meaning the tip can be deformed in order to enter gently through a clear corneal incision, after which it assumes its original round shape within the anterior chamber. **Complete, round capsulotomies are accomplished in a few milliseconds.**


## HOW DOES IT WORK?

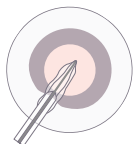


**ONE:** Capsulotomy tip extends into elongated shape, surgeon extends pushrod for entry through corneal incision.

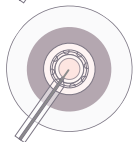

**TWO:** Surgeon retracts the pushrod, which allows the capsulotomy tip to naturally return to a circular shape within the anterior chamber. The surgeon then centers the suction cup either on the pupil or the visual axis of the patient.

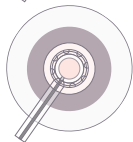


**THREE:** After the desired capsulotomy placement, a small amount of suction is applied through the silicone suction cup to appose the nitinol capsulotomy ring to the lens capsule.

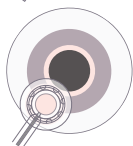

**FOUR:** At the surgeon's command, the capsulotomy is performed in milliseconds. Suction is relieved as the surgeon slides the sliding mechanism to the back position and ZEPTO™ is retrieved through a corneal incision.

# YOUR PATIENTS

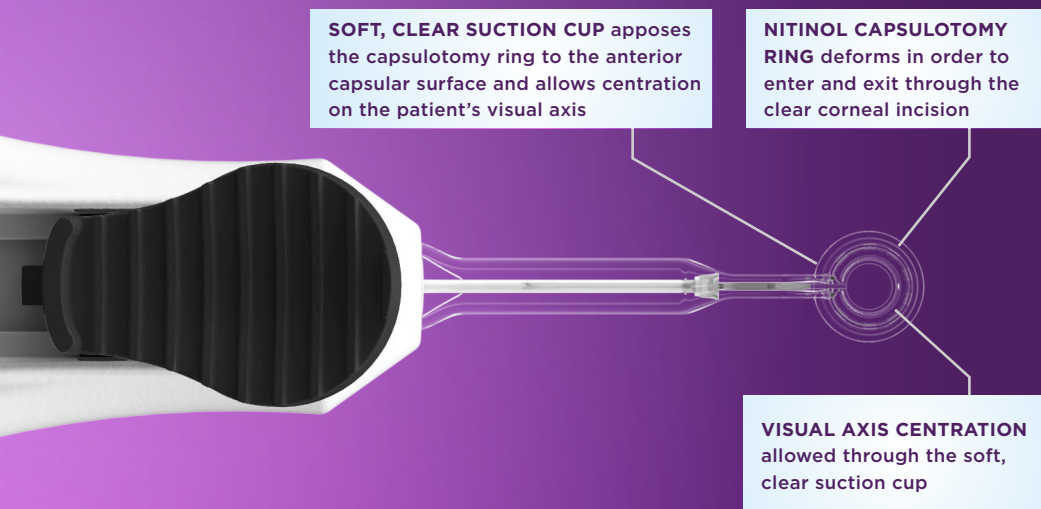

## SAFE AND EFFECTIVE

FDA Clinical Trial completed enrollment in December 2016.

Safety testing in rabbit eyes has included temperature measurements using thermocouples in the eye during ZEPTO™ capsulotomy. Only very brief transient temperature changes of 1° to 2°C lasting several seconds were noted adjacent to the ZEPTO™ tip and the corneal endothelium.<sup>1</sup>

Zonular safety has also been investigated using Miyake-Apple imaging during ZEPTO™ capsulotomy in rabbit eyes. Very little stress on the zonules was observed during application of suction, capsulotomy, and suction release.<sup>1</sup>

A complete assessment of safety and performance has been performed in live rabbits, which showed no corneal endothelial cell damage, inflammatory response, or gross and microscopic tissue damage due to ZEPTO™ use.<sup>1</sup>

## UNIQUE FEATURES

ZEPTO™ capsulotomies have a **capsular edge tear strength 2 to 4 times greater** than that of continuous curvilinear capsulorhexis and femtosecond laser capsulotomies.<sup>2</sup>

ZEPTO™ allows **intraoperative capsulotomy centration on the patient's visual axis** for an optimized outcome.

**Mynosys**

# A NOVEL DISPOSABLE CAPSULOTOMY DEVICE

ZEPTO™ allows you to create precise, visually centered, and stronger capsulotomies in milliseconds.

## WHAT IS ZEPTO™?

- » Instantaneous capsulotomy
- » Precision centration on the visual axis
- » Automated capsulotomy technology = ease of use
- » Works with small pupils
- » Not affected by corneal scars or irregularities
- » Easy practice adoption
- » Seamless integration into surgical sequence
- » Ideal for difficult cases like weak zonules and dense cataracts

For more information, visit [www.mynosys.com](http://www.mynosys.com)



**REFERENCES:** **1.** Chang DF, et al. Precision pulse capsulotomy – preclinical safety and performance of a new capsulotomy technology. *Ophthalmology*. In press. **2.** Thompson V, et al. Comparison of manual, femtosecond laser, and precision pulse capsulotomy edge tear strength in paired human cadaver eyes. *Ophthalmology*. In press.

ZEPTO™ is a trademark of Mynosys Cellular Devices, Inc. • 12763 Rev B

**Mynosys**